you must go on

- EN

- Liposuction Corsets are recommended for post-operative period after liposuction operations such as on abdominal, hip, back, thigh, calf, arm etc. Helps to reduce swelling and scar, to shape the body contours and also to speed up the healing process after liposuction operations. Provides controlled pressure to the area where it is applied.
 - During the use of the corset after liposuction surgery, ancillary products (such as pads, boards) are also given in case extra support is needed to ensure that the corset fits the body contours better. These ancillary supports are used together with the relevant liposuction corset. Different ancillary supports which are used on abdomen, back, waist or in areas where uniform healing is required, help the corset to do function better after liposuction surgery. It also acts as a soft barrier between the corset and the skin which is sensitive after operation.
- TR

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EN Storage

Must be kept in the closed packaging and without being exposed to sun, dust or humidity.

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EN FEATURES

The product is a foam compression support with adhesive tape. Specially designed to be used under corset after liposuction surgeries. Easy to apply under the corset and comfortable to use. If necessary, it can be cut and used in shapes and sizes appropriate to the area where it will be applied.

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Dimensions (-): 200mm * 300mm * 10mm

EN WASHING/CLEANING

The product cannot be cleaned and reused. Therefore, if the need continues, the product should be replaced with a new one rather than cleaning and reusing it.

TR YIKAMA/TEMİZLEME

The product cannot be cleaned and reused. Therefore, if the need continues, the product should be replaced with a new one rather than cleaning and reusing it.

FR LAVAGE/NETTOYAGE

The product cannot be cleaned and reused. Therefore, if the need continues, the product should be replaced with a new one rather than cleaning and reusing it.

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EN HOW TO APPLY

EN As seen in the drawing, it is applied under the corset in the needed area. It is applied with the adhesive part on the skin. This feature helps the product fit snugly into the area where it is applied and remain stable. If necessary, it can be cut and can be used in the shapes and sizes appropriate to the area where it will be applied.

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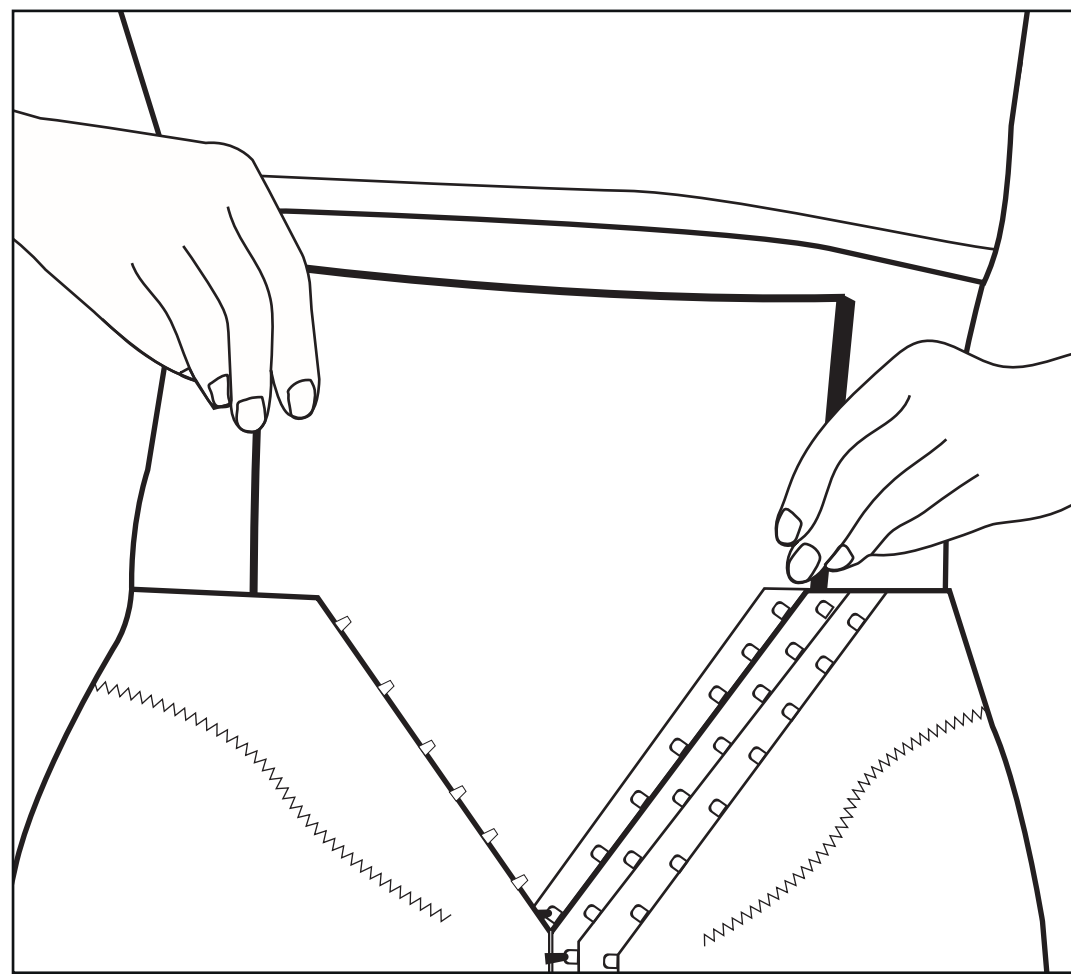
EN WARNING

It is used as per doctor's advice. After the surgery, the use of corset and ancillary supports together with the corset should be done under regular doctor control and supervision. It is not used on open wounds. In case of skin rash or unexpected cases, consult your doctor. Stop using the product, if you have sensitivity or allergy against the materials which the product contains. Do not contact with eyes. If eye irritation occurs due to contact, wash the entire eye with plenty of water and consult a doctor.

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